

anti- ABLIM1 antibody

Product Information

Catalog No.: FNab00054

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

May act as scaffold protein(By similarity). May play a role in the development of the retina. Has been suggested to play a role in axon guidance.

Immunogen information

Immunogen: actin binding LIM protein 1

Synonyms: ABLIM, abLIM 1, ABLIM1, actin binding LIM protein 1,

DKFZp781D0148, FLJ14564, KIAA0059, LIMAB1, LIMATIN

Observed MW: 88 kDa UniprotID: 014639

Application

Reactivity: Human, Mouse, Rat

Tested Application: ELISA, WB, IHC

Recommended dilution: WB: 1:500-1:5000; IHC: 1:20-1:200

Image:


1

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 www.fn-test.com



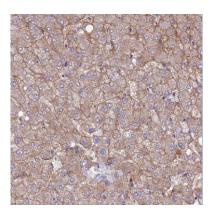

Immunohistochemistry of paraffin-embedded human liver using FNab00054(ABLIM1 antibody) at dilution of 1:100

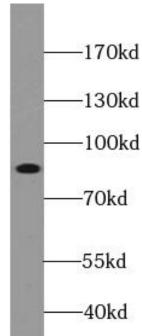

HEK-293 cells were subjected to SDS PAGE followed by western blot with FNab00054(ABLIM1 antibody) at dilution of 1:1000