

anti- AIRE antibody

Product Information

Catalog No.: FNab00241

Size: 100µg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

Transcriptional regulator that binds to DNA as a dimer or as a tetramer, but not as a monomer. Binds to G-doublets in an A/T-rich environment; the preferred motif is a tandem repeat of 5'-. ATTGGTTA-3' combined with a 5'-TTATTA-3' box. Binds to nucleosomes(By similarity). Binds to chromatin and interacts selectively with histone H3 that is not methylated at 'Lys-4', not phosphorylated at 'Thr-3' and not methylated at 'Arg-2'. Functions as a sensor of histone H3 modifications that are important for the epigenetic regulation of gene expression. Functions as a transcriptional activator and promotes the expression of otherwise tissue-specific self-antigens in the thymus, which is important for self tolerance and the avoidance of autoimmune reactions.

Immunogen information

Immunogen: autoimmune regulator

Synonyms: AIRE, AIRE1, APECED, APECED protein, APS1, APSI, autoimmune

regulator, PGA1

Observed MW: 60 kDa UniprotID: 043918

Application

1

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 www.fn-test.com




Reactivity: Human, Mouse


Tested Application: ELISA, WB, IHC

Recommended dilution: WB: 1:500-1:5000; IHC: 1:20-200

Image:

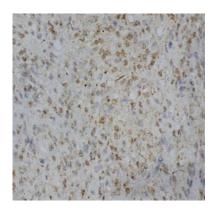

Immunohistochemistry of paraffin-embedded human stomach slide using FNab00241(AIRE Antibody) at dilution of 1:50

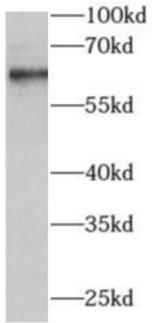

L02 cells were subjected to SDS PAGE followed by western blot with FNab00241(AIRE Antibody) at dilution of 1:1000

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel: (0086)027-87384275 Fax: (0086)027-87800889 www.fn-test.com