

anti- ALG6 antibody

Product Information

Catalog No.:	FNab00308
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Adds the first glucose residue to the lipid-linked oligosaccharide precursor for N-linked glycosylation. Transfers glucose from dolichyl phosphate glucose(Dol-P-Glc) onto the lipid-linked oligosaccharide Man(9)GlcNAc(2)-PP-Dol.

Immunogen information

Immunogen:	asparagine-linked glycosylation 6, alpha-1,3-glucosyltransferase homolog(<i>S. cerevisiae</i>)
Synonyms:	ALG6
Observed MW:	58 kDa
UniprotID :	Q9Y672

Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB
Recommended dilution:WB:	1:200-1:2000


Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

Image:

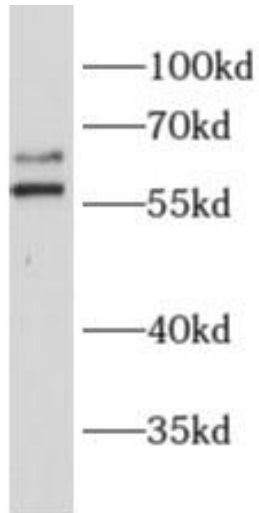

HeLa cells were subjected to SDS PAGE followed by western blot with FNab00308(ALG6 antibody) at dilution of 1:600

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com