

anti- Angiotensinogen antibody

Product Information

Catalog No.:	FNab00393
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	\geq 95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Angiotensinogen is a precursor of angiotensin II(Ang II), is expressed and synthesized largely in the liver and is cleaved by the enzyme renin in response to lowered blood pressure. It has a key role in mediating vascular constriction and regulating salt and fluid homeostasis. The resulting product, angiotensin I, is then cleaved by angiotensin converting enzyme(ACE) to generate the physiologically active enzyme angiotensin II. Mutations in this gene are associated with susceptibility to essential hypertension, and can cause renal tubular dysgenesis, a severe disorder of renal tubular development. Defects in this gene have also been associated with non-familial structural atrial fibrillation, and inflammatory bowel disease.

Immunogen information

Immunogen:	angiotensinogen(serpin peptidase inhibitor, clade A, member 8)
Synonyms:	AGT, Ang I Angiotensin 2, Ang II Angiotensin 3, Ang III, Angiotensin I, Angiotensin II, Angiotensin III, Angiotensinogen, ANHU, Des Asp[1] angiotensin II, FLJ92595, FLJ97926, Serpin A8, SERPINA8
Observed MW:	53 kDa
UniprotID :	P01019

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>



Application

Reactivity:Human, MouseTested Application:ELISA, WB, IHC, FCRecommended dilution:WB: 1:200-1:2000; IHC: 1:20-1:200Image:

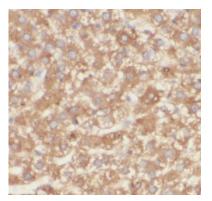

Immunohistochemistry of paraffin-embedded human liver slide using FNab00393(AGT Antibody) at dilution of 1:50

HepG2 cells were subjected to SDS PAGE followed by western blot with FNab00393(AGT antibody) at dilution of 1:600

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>