

# anti-Annexin VIII antibody

### **Product Information**

Catalog No.:	FNab00440
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	$\geq$ 95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)
De de successi d	

#### Background

None

## **Immunogen information**

Immunogen:	annexin A8-like 2
Synonyms:	Annexin 8, Annexin A8, annexin A8 like 2, Annexin VIII, ANX8, ANXA8, ANXA8L2, bA145E20.2, FLJ32754, FLJ54151, VAC beta, Vascular anticoagulant beta
Observed MW:	37 kDa
UniprotID :	P13928

## Application


Reactivity:	Human, Mouse, Rat	
Tested Application:	ELISA, WB, IHC	
Recommended dilution: WB: 1:200-1:1000; IHC: 1:50-1:500		
Image:		

#### Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206) Tel :( 0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>

1



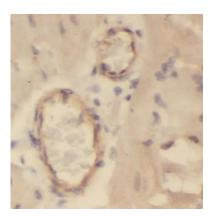

Immunohistochemistry of paraffin-embedded human skin tissue slide using FNab00440( Annexin VIII antibody) at dilution of 1:200

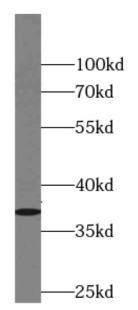

HepG2 cells were subjected to SDS PAGE followed by western blot with FNab00440(ANXA8 antibody) at dilution of 1:100

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :( 0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>