

anti- ATG16L1 antibody

Product Information

Catalog No.: FNab00669

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

The protein encoded by this gene is part of a large protein complex that is necessary for autophagy, the major process by which intracellular components are targeted to lysosomes for degradation. Defects in this gene are a cause of susceptibility to inflammatory bowel disease type 10 (IBD10). Several transcript variants encoding different isoforms have been found for this gene.

Immunogen information

Immunogen: ATG16 autophagy related 16-like 1

Synonyms: APG16L
Observed MW: 68 kDa
Uniprot ID: Q676U5

Application

Reactivity: Human, Mouse

Tested Application: ELISA, WB, IF, IHC

Recommended dilution: WB: 1:500 - 1:2000; IHC: 1:50 - 1:200; IF: 1:50 - 1:200

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>


1



Image:

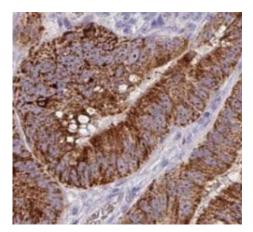

Immunohistochemistry of paraffin-embedded human colon using FNab00669(ATG16L1 antibody) at dilution of 1:100

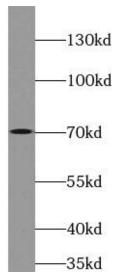

MCF7 cells were subjected to SDS PAGE followed by western blot with FNab00669(ATG16L1 antibody) at dilution of 1:1000