

anti- BLNK antibody

Product Information

Catalog No.: FNab00908

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

Functions as a central linker protein, downstream of the B-cell receptor(BCR), bridging the SYK kinase to a multitude of signaling pathways and regulating biological outcomes of B-cell function and development. Plays a role in the activation of ERK/EPHB2, MAP kinase p38 and JNK.Modulates AP1 activation.Important for the activation of NF-kappa-B and NFAT. Plays an important role in BCR-mediated PLCG1 and PLCG2 activation and Ca(2+) mobilization and is required for trafficking of the BCR to late endosomes. However, does not seem to be required for pre-BCR-mediated activation of MAP kinase and phosphatidyl-inositol 3(PI3) kinase signaling. May be required for the RAC1-JNK pathway.Plays a critical role in orchestrating the pro-B cell to pre-B cell transition.May play an important role in BCR-induced B-cell apoptosis.

Immunogen information

Immunogen: B-cell linker
Synonyms: BASH, SLP65
Observed MW: 68 kDa, 70 kDa

UniprotID: Q8WV28

Application

1

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 www.fn-test.com




Reactivity: Human, Mouse, Rat


Tested Application: ELISA, WB, IHC

Recommended dilution: WB: 1:500-1:2000; IHC: 1:20-1:200

Image:

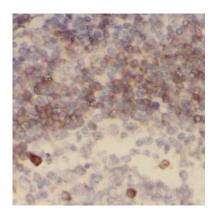

Immunohistochemistry of paraffin-embedded human tonsillitis tissue slide using FNab00908(BLNK Antibody) at dilution of 1:400 heat mediated antigen retrieved with Tris-EDTA buffer(pH9).

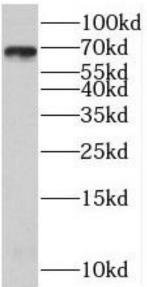

Raji cells were subjected to SDS PAGE followed by western blot with FNab00908(BLNK Antibody) at dilution of 1:1000