

anti- BZW1 antibody

Product Information

Catalog No.:	FNab00995
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Enhances histone H4 gene transcription but does not seem to bind DNA directly.

Immunogen information

Immunogen:	basic leucine zipper and W2 domains 1
Synonyms:	BZAP45, KIAA0005
Observed MW:	48kd
UniprotID :	Q7L1Q6

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC, IF
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200; IF: 1:20-1:200
Image:	

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

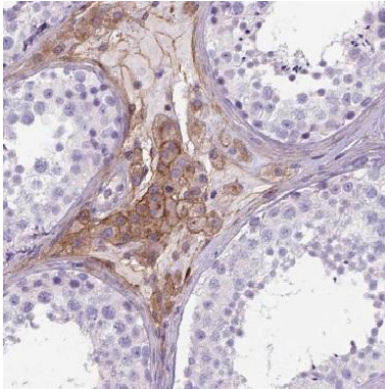

Immunohistochemistry of paraffin-embedded human testis using FNab00995(BZW1 antibody) at dilution of 1:100

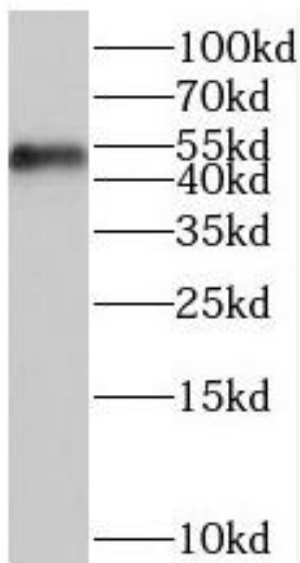

MCF7 cells were subjected to SDS PAGE followed by western blot with FNab00995(BZW1 antibody) at dilution of 1:500