

anti- CD63 antibody

Product Information

Catalog No.: FNab01490

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

The protein encoded by this gene is a member of the transmembrane 4 superfamily, also known as the tetraspanin family. Most of these members are cell-surface proteins that are characterized by the presence of four hydrophobic domains. The proteins mediate signal transduction events that play a role in the regulation of cell development, activation, growth and motility. The encoded protein is a cell surface glycoprotein that is known to complex with integrins. It may function as a blood platelet activation marker. Deficiency of this protein is associated with Hermansky-Pudlak syndrome. Also this gene has been associated with tumor progression. Alternative splicing results in multiple transcript variants encoding different protein isoforms.

Immunogen information

Immunogen: CD63 molecule

Synonyms: MLA1, TSPAN30

Observed MW: 26-60 kDa
Uniprot ID: P08962

Application

Reactivity: Human, Mouse, Rat

1

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)


Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>




Tested Application: ELISA, IHC, WB

Recommended dilution: WB: 1:500 - 1:2000; IHC: 1:50 - 1:200

Image:

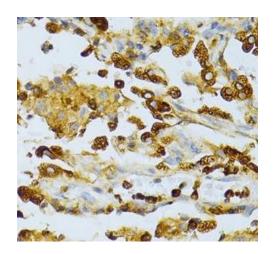

Immunohistochemistry of paraffin-embedded human lung cancer using FNab01490(CD63 Antibody) at dilution of 1:100 heat mediated antigen retrieved with Tris-EDTA buffer(pH9).

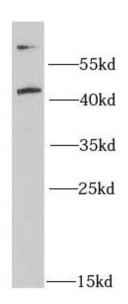

Mouse spleen were subjected to SDS PAGE followed by western blot with FNab01490(CD63 Antibody) at dilution of 1:1000