

anti- COX4I2 antibody

Product Information

Catalog No.:	FNab01899
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

This protein is one of the nuclear-coded polypeptide chains of cytochrome c oxidase, the terminal oxidase in mitochondrial electron transport.

Immunogen information

Immunogen:	cytochrome c oxidase subunit IV isoform 2(lung)
Synonyms:	COX4L2
Observed MW:	17-18 kDa
UniprotID :	Q96KJ9

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC, IF, IP
Recommended dilution:	WB: 1:500-1:2000; IP: 1:200-1:1000; IHC: 1:20-1:200; IF: 1:20-1:200
Image:	

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275


Fax: (0086)027-87800889 www.fn-test.com

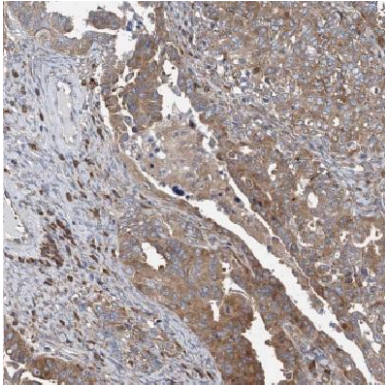

Immunohistochemistry of paraffin-embedded human ovary tumor using FNab01899(COX4I2 antibody) at dilution of 1:50

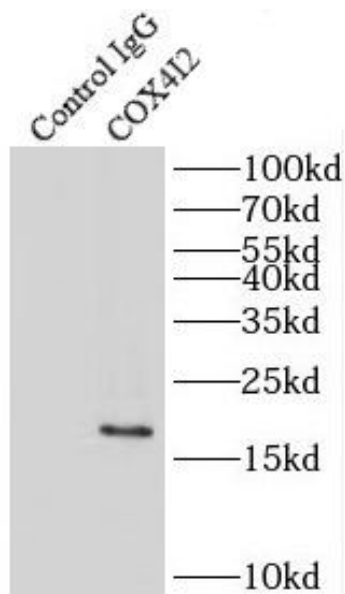

IP Result of anti-COX4I2 (IP:FNab01899, 3ug; Detection:FNab01899 1:1000) with HepG2 cells lysate 8000ug.

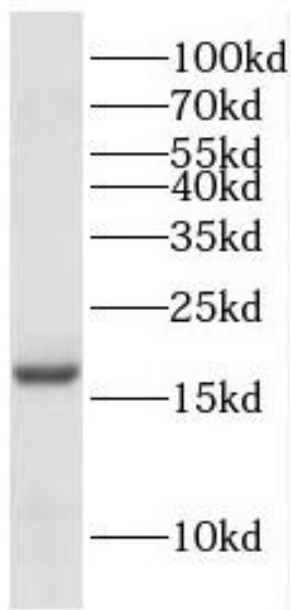

MCF7 cells were subjected to SDS PAGE followed by western blot with FNab01899(COX4I2 antibody) at dilution of 1:1000