

anti- NCOA5 antibody

Product Information

Catalog No.:	FNab05590
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Nuclear receptor coregulator that can have both coactivator and corepressor functions. Interacts with nuclear receptors for steroids(ESR1 and ESR2) independently of the steroid binding domain(AF-2) of the ESR receptors, and with the orphan nuclear receptor NR1D2. Involved in the coactivation of nuclear steroid receptors(ER) as well as the corepression of MYC in response to 17-beta-estradiol(E2).

Immunogen information

Immunogen:	nuclear receptor coactivator 5
Synonyms:	KIAA1637
Observed MW:	66 kDa
Uniprot ID :	Q9HCD5

Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IF
Recommended dilution:	WB: 1:500-1:2000; IF: 1:20-1:200


Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

Image:

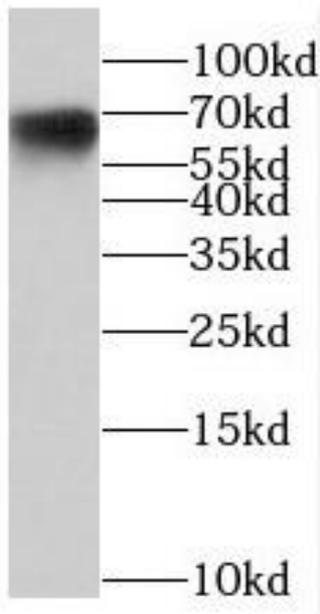

HEK-293 cells were subjected to SDS PAGE followed by western blot with FNab05590(NCOA5 antibody) at dilution of 1:600

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com