

anti- NOMO1 antibody

Product Information

Catalog No.:	FNab05791
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

May antagonize Nodal signaling.

Immunogen information

Immunogen:	NODAL modulator 1
Synonyms:	PM5
Observed MW:	135 kDa
UniprotID :	Q15155

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IF, IHC
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200; IF: 1:20-1:200
Image:	

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

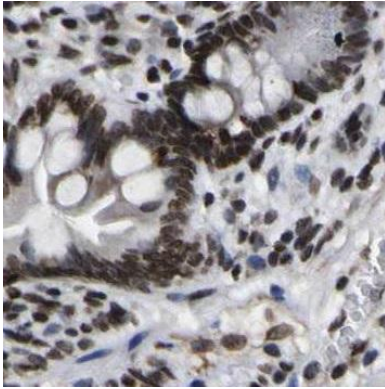

Immunohistochemistry of paraffin-embedded human colon slide using FNab05791(NOMO1 Antibody) at dilution of 1:50

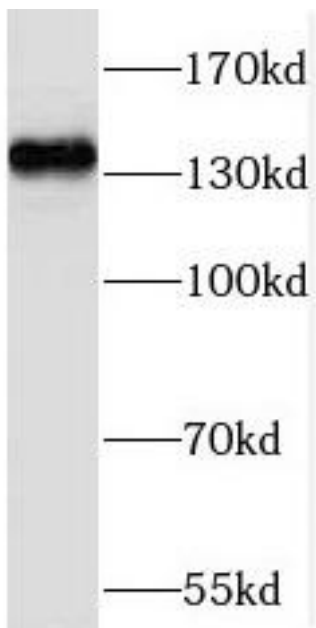

mouse colon tissue were subjected to SDS PAGE followed by western blot with FNab05791(NOMO1 antibody) at dilution of 1:500