

anti- OLFM1 antibody

Product Information

Catalog No.:	FNab05983
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Contributes to the regulation of axonal growth in the embryonic and adult central nervous system by inhibiting interactions between RTN4R and LINGO1. Inhibits RTN4R-mediated axon growth cone collapse(By similarity). May play an important role in regulating the production of neural crest cells by the neural tube(By similarity). May be required for normal responses to olfactory stimuli(By similarity).

Immunogen information

Immunogen:	olfactomedin 1
Synonyms:	NOE1, NOEL1
Observed MW:	None
Uniprot ID :	Q99784

Application

Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, IHC
Recommended dilution:IHC:	1:20-1:200

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889

www.fn-test.com

Image:

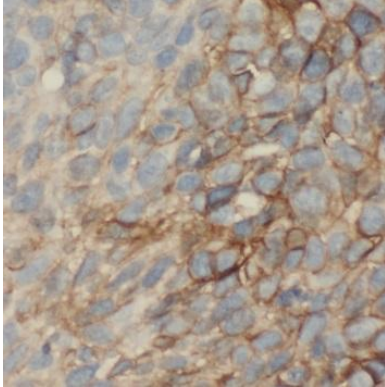

Immunohistochemistry of paraffin-embedded human lung cancer tissue using FNab05983(OLFM1 antibody) at dilution of 1:50