

anti- RRBP1 antibody

Product Information

Catalog No.:	FNab07491
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Acts as a ribosome receptor and mediates interaction between the ribosome and the endoplasmic reticulum membrane.

Immunogen information

Immunogen:	ribosome binding protein 1 homolog 180kDa(dog)
Synonyms:	KIAA1398
Observed MW:	180 kDa, 109 kDa
UniprotID :	Q9P2E9

Application


Reactivity:	Human, Mouse
Tested Application:	ELISA, WB, IHC, IF
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200; IF: 1:20-1:200
Image:	

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

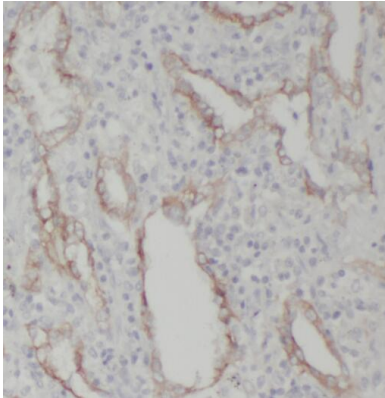

Immunohistochemistry of paraffin-embedded human lung cancer tissue slide using FNab07491(RRBP1 Antibody) at dilution of 1:50

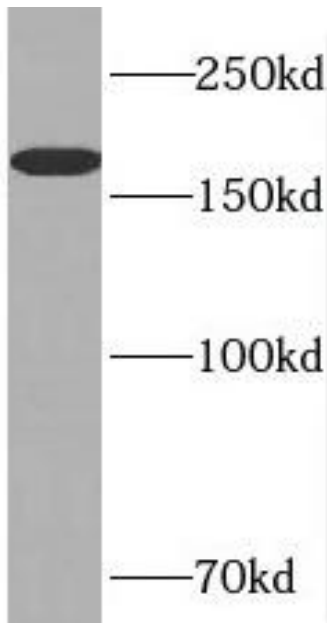

mouse pancreas tissue were subjected to SDS PAGE followed by western blot with FNab07491(RRBP1 antibody) at dilution of 1:800