

anti- TMEFF2 antibody

Product Information

Catalog No.:	FNab08742
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

May be a survival factor for hippocampal and mesencephalic neurons. The shedded form up-regulates cancer cell proliferation, probably by promoting ERK1/2 phosphorylation.

Immunogen information

Immunogen:	transmembrane protein with EGF-like and two follistatin-like domains 2
Synonyms:	HPP1, TENB2, TPEF
Observed MW:	45-48 kDa
UniprotID :	Q9UIK5

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200
Image:	

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

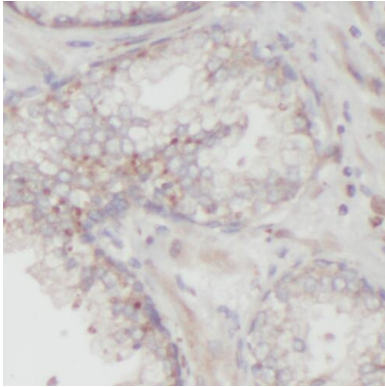

Immunohistochemistry of paraffin-embedded human prostate cancer tissue slide using FNab08742(TMEFF2 Antibody) at dilution of 1:100

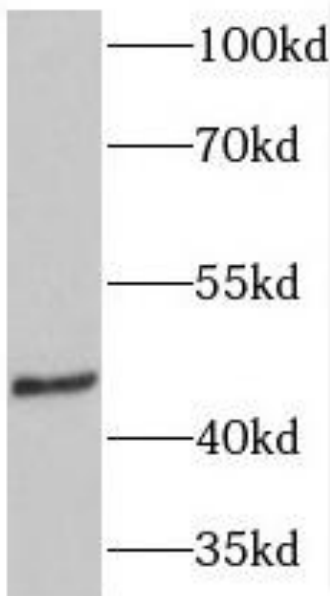

SH-SY5Y cells were subjected to SDS PAGE followed by western blot with FNab08742(TMEFF2 antibody) at dilution of 1:300