

anti- CAV3 antibody

Product Information

Catalog No.: FNab10145

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

This gene encodes a caveolin family member, which functions as a component of the caveolae plasma membranes found in most cell types. Caveolin proteins are proposed to be scaffolding proteins for organizing and concentrating certain caveolin-interacting molecules. Mutations identified in this gene lead to interference with protein oligomerization or intra-cellular routing, disrupting caveolae formation and resulting in Limb-Girdle muscular dystrophy type-1C (LGMD-1C), hyperCKemia or rippling muscle disease (RMD). Alternative splicing has been identified for this locus, with inclusion or exclusion of a differentially spliced intron. In addition, transcripts utilize multiple polyA sites and contain two potential translation initiation sites.

Immunogen information

Immunogen: caveolin 3

Synonyms: CAV3, caveolin 3, Caveolin-3, LGMD1C, LQT9, M caveolin, VIP 21,

VIP21

Observed MW: 17kDa Uniprot ID: P56539

Application

1

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>




Reactivity: Human, Mouse


Tested Application: ELISA, WB, IHC

Recommended dilution: WB: 1:500 - 1:2000; IHC: 1:50 - 1:200

Image:

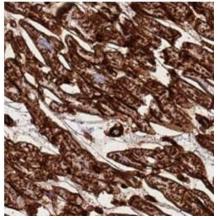

Immunohistochemistry of paraffin-embedded human heart using FNab10145(CAV3 antibody) at dilution of 1:100

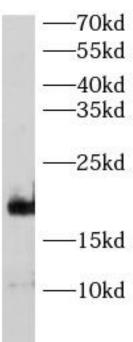

mouse heart were subjected to SDS PAGE followed by western blot with FNab10145(CAV3 antibody) at dilution of 1:1000