

anti- HYOU1 antibody

Product Information

Catalog No.:	FNab10247
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	\geq 95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

The protein encoded by this gene belongs to the heat shock protein 70 family. This gene uses alternative transcription start sites. A cis-acting segment found in the 5' UTR is involved in stress-dependent induction, resulting in the accumulation of this protein in the endoplasmic reticulum (ER) under hypoxic conditions. The protein encoded by this gene is thought to play an important role in protein folding and secretion in the ER. Since suppression of the protein is associated with accelerated apoptosis, it is also suggested to have an important cytoprotective role in hypoxia-induced cellular perturbation. This protein has been shown to be up-regulated in tumors, especially in breast tumors, and thus it is associated with tumor invasiveness. This gene also has an alternative translation initiation site, resulting in a protein that lacks the N-terminal signal peptide. This signal peptide-lacking protein, which is only 3 amino acids shorter than the mature protein in the ER by a carboxy-terminal peptide sequence and to mitochondria by an amino-terminal targeting signal. Alternative splicing results in multiple transcript variants.

Immunogen information

Immunogen:	Hypoxia up-regulated protein 1
Synonyms:	GRP-170, Grp170, HSP12A, ORP-150, ORP150
Observed MW:	110kDa, 150kDa

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889


www.fn-test.com




Uniprot ID : Q9Y4L1

Application

Reactivity:Human, MouseTested Application:ELISA, WB, IHCRecommended dilution:WB: 1:500 - 1:2000; IHC: 1:50 - 1:200Image:

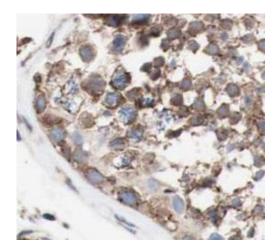

Immunohistochemistry of paraffin-embedded mouse testis tissue slide using FNab10247(HYOU1 Antibody) at dilution of 1:100

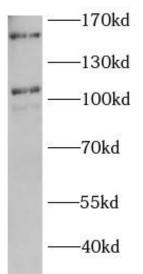

MCF7 cells were subjected to SDS PAGE followed by western blot with FNab10247(HYOU1 antibody) at dilution of 1:1000

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889

www.fn-test.com