


Recombinant Human MAP3K5

Catalogue No.:	P2610
Species:	Human
Uniprot ID:	Q99683
Expression Region:	1180-1374
Host:	E.Coli
Tags:	N-terminal His-IF2DI Tag
Molecular Weight:	42.9 kDa under reducing conditions
Purity:	Greater than 90% as determined by SDS-PAGE
Formulation:	Lyophilized from a 0.2 μ m filtered solution in 10 mM Hepes, 500 mM NaCl with 5% trehalose, pH7.4
Reconstitution:	Centrifuge the vial prior to opening, reconstitute in sterile distilled water or buffer of interest to a concentration of 0.1-1mg/ml by gently pipetting 2-3 times, don't vortex.
Storage:	The lyophilized protein is stable at -20°C for up to 1 year, the protein solution can be stored at 2-8°C for up to 1 week. For extended storage, it is recommended to further dilute in a buffer containing 0.1% BSA (make sure BSA will NOT interfere with your experiment) and store in working aliquots at -20°C to -80°C for 6 months. Avoid repeated freeze/thaw cycle.

Synonyms: Apoptosis signal regulating kinase 1, Apoptosis signal-regulating kinase 1, ASK 1, ASK-1, ASK1, M3K5, M3K5_HUMAN, MAP/ERK kinase kinase 5, MAP3K5, MAPK/ERK kinase kinase 5, MAPKKK5, MEK kinase 5, MEKK 5, MEKK5, Mitogen activated protein kinase kinase kinase 5, Mitogen-activated protein kinase kinase kinase 5

SDS-PAGE:

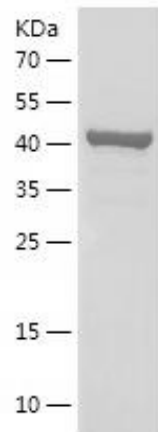

Safety Note: This product is intended for research and manufacturing uses only. It is not a diagnostic device. Product degradation will result from multiple freeze/thaw cycles. It is suggested that the antigen be stored in use size aliquots and thawed just prior to use. This material has been inactivated, however as with all biological materials, it should be handled as potentially infectious. The user assumes all responsibility for care, custody and control of the material, including its disposal, in accordance with all regulations.