


Recombinant Mouse IL-22

Catalogue No.:	P6731
Species:	Mouse
Uniprot ID:	Q9JJY9
Expression Region:	34-179
Host:	E.Coli
Tags:	N-terminal His Tag
Molecular Weight:	18.5 kDa under reducing conditions
Purity:	Greater than 95% as determined by SDS-PAGE
Formulation:	Lyophilized from a 0.2 μ m filtered solution in 10 mM Hepes, 500 mM NaCl with 5% trehalose, pH7.4
Reconstitution:	Centrifuge the vial prior to opening, reconstitute in sterile distilled water or buffer of interest to a concentration of 0.1-1mg/ml by gently pipetting 2-3 times, don't vortex.
Storage:	The lyophilized protein is stable at -20°C for up to 1 year, the protein solution can be stored at 2-8°C for up to 1 week. For extended storage, it is recommended to further dilute in a buffer containing 0.1% BSA (make sure BSA will NOT interfere with your experiment) and store in working aliquots at -20°C to -80°C for 6 months. Avoid repeated freeze/thaw cycle.

Synonyms: Cytokine Zcyto18, IL 10 related T cell derived inducible factor, IL 21, IL 22, IL D110, IL TIF, IL-10-related T-cell-derived-inducible factor, IL-22, IL-TIF, IL21, IL22, IL22_HUMAN, ILD110, ILTIF, Interleukin 10 related T cell derived inducible factor

SDS-PAGE:

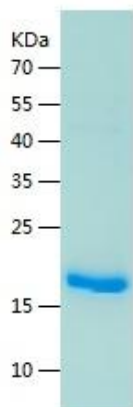

Safety Note: This product is intended for research and manufacturing uses only. It is not a diagnostic device. Product degradation will result from multiple freeze/thaw cycles. It is suggested that the antigen be stored in use size aliquots and thawed just prior to use. This material has been inactivated, however as with all biological materials, it should be handled as potentially infectious. The user assumes all responsibility for care, custody and control of the material, including its disposal, in accordance with all regulations.