


Recombinant Human SAE1

Catalogue No.:	P6890
Species:	Human
Uniprot ID:	Q9UBE0
Expression Region:	1-346
Host:	E.Coli
Tags:	N-terminal His Tag
Molecular Weight:	38 kDa under reducing conditions
Purity:	Greater than 95% as determined by SDS-PAGE
Formulation:	Lyophilized from a 0.2 μ m filtered solution in 10 mM Hepes, 500 mM NaCl with 5% trehalose, pH7.4.
Reconstitution:	Centrifuge the vial prior to opening, reconstitute in sterile distilled water or buffer of interest to a concentration of 0.1-1mg/ml by gently pipetting 2-3 times, don't vortex.
Storage:	The lyophilized protein is stable at -20°C for up to 1 year, the protein solution can be stored at 2-8°C for up to 1 week. For extended storage, it is recommended to further dilute in a buffer containing 0.1% BSA (make sure BSA will NOT interfere with your experiment) and store in working aliquots at -20°C to -80°C for 6 months. Avoid repeated freeze/thaw cycle.

Synonyms: Activator of SUMO1, AOS1, HSPC140, Sae1, SAE1_HUMAN, Sentrin/SUMO activating protein AOS1, SUA1, SUMO 1 activating enzyme E1 N subunit, SUMO 1 activating enzyme subunit 1, SUMO-activating enzyme subunit 1, Ubiquitin like protein SUMO1 activating enzyme, Ubiquitin-like 1-activating enzyme E1A, UBL E1A, UBLE1A

Image:

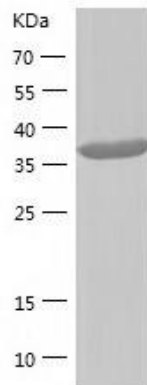

Safety note: This product is intended for research and manufacturing uses only. It is not a diagnostic device. Product degradation will result from multiple freeze/thaw cycles. It is suggested that the antigen be stored in use size aliquots and thawed just prior to use. This material has been inactivated, however as with all biological materials, it should be handled as potentially infectious. The user assumes all responsibility for care, custody and control of the material, including its disposal, in accordance with all regulations.